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OPERATING INSTRUCTIONS

Venen Engel 3-in-1 Glow Brush

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Venen Engel 3-in-1 Glow Brush

PRELIMINARY REMARK

Congratulations on purchasing this device! You will enjoy using your Venen Engel Product! If you are satisfied with your Venen Engel, we would be delighted if you would recommend it to others.

This manual contains details that users should observe in order to prevent hazards and use the device safely.

If you have any questions about use, accessories, or suggestions, please do not hesitate to contact us!

Mail service@venenengel.de

Web www.venenengel.de

BEFORE YOU START

- Please read the operating instructions carefully before using the device.
- Keep this instruction manual in a safe place so that you can refer to it when necessary
- The Venen Engel product must not be used in certain health conditions. Please read the list of contraindications carefully.

IMPORTANT SAFETY INSTRUCTIONS AND WARNINGS

It is important that you read all warnings and precautions in this manual. They are for your safety, to prevent injury, and to avoid damage to the device.

- Please read these instructions for use carefully before using the device.
- Keep this manual for future reference.
- Keep the device and its accessories out of reach of children.
- The device is intended for capable adults only.
- Do not use the device if it is damaged.
- Do not use the device if the charging cable is damaged.
- Do not open or disassemble the device.
- Only charge the device using the supplied USB cable.
- Do not look directly into the LEDs during operation or shine the light directly into the eyes of other persons.
- The device has been tested for resistance to **temporary water exposure**. However, to ensure safe operation and product longevity, the main unit must not be exposed to direct water jets or continuous running water.

- If the device is used in a shower environment, avoid direct contact with the water stream and keep the main unit away from excessive moisture.
- Do not immerse the main unit in water.
- The massage heads may be removed from the main unit for cleaning. Clean them separately with water if necessary and allow them to dry completely before reattaching them.
- Do not use the device near hair to avoid entanglement.
- Do not operate the device continuously for more than 30 minutes.
- Do not use the device on the same body area for more than 15 minutes.
- Turn off the device before changing the massage heads.

SCOPE OF DELIVERY

Body Massager


USB charging cable

Roller massage head

Soft brush
massage head

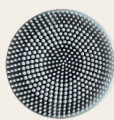Brush massage head

TECHNICAL DATA

Modell	PL-674
Rated voltag	5V DC
Power	16 W
Recommended operating time	15 minutes per body area. No longer than 30 minutes of continuous use
Power supply	Rechargeable battery (USB charging)
Packaging Dimensions	14,5 x 2,8 x 9,6cm
Weight	0,5kg

EXPLANATION OF SYMBOLS

Manufacturer information

Conformity with European standards

Made in China

Electrical product that contains a battery. Please follow disposal instructions.

Read manual before use

IPX7

Product is protected against **temporary submersion**, not continuous, high-pressure, or deep water.

Serial number

Inner and outer packaging made of mixed paper

USB port

Disposal instructions France for electronic products

Disposal instructions France for the packaging

PRODUCT DESCRIPTION

This compact, rechargeable body massager is designed to complement a daily wellness and self-care routine. This versatile device delivers a gentle, relaxing surface massage for both muscles and skin. It provides a comforting experience to help muscles relax after exercise or daily physical activity across targeted areas of the body, including the back, shoulders, arms, legs, and feet.

To customize the massage experience, the 3-in-1 Glow Brush features interchangeable, rotating massage heads that deliver a variety of massage styles and sensations. Certain operating modes also feature integrated red LED lights that illuminate during use to enhance the application experience. In addition to muscle relaxation, the specialized attachments may be integrated into a cosmetic skincare routine for gentle body cleansing and exfoliation.

The Venen Engel 3-in-1 Glow Brush is intended strictly for cosmetic, personal care, and wellness purposes. It is not a medical device and is not intended to diagnose, treat, cure, or prevent any medical condition.

INTENDED USE

The **Venen Engel 3-in-1 Glow Brush** is a rechargeable device intended for cosmetic and wellness massage of the body.

The device is designed to provide a relaxing surface massage of muscles and skin using rotating massage heads. It may help users relax muscles after physical activity or daily strain. The device is intended to be used on external body areas, including:

- back
- shoulders
- arms
- legs
- feet

The interchangeable massage attachments allow different massage sensations and may also be used for cosmetic skin care purposes, such as gentle body cleansing or exfoliation.

Target user

The product is intended for adult users.

CONTRAINDICATIONS

PLEASE CONSULT A PHYSICIAN BEFORE USING THE DEVICE IF YOU:

- are pregnant
- have a serious heart condition
- have injuries or skin irritations in the treatment areas.

If you are unsure if the device is suitable for you, consult a physician before use.

Stop using the device immediately if you feel discomfort or pain.

IMPROPER USE

Any use other than the intended use is considered improper.

Examples of foreseeable misuse include:

- Use on the head or near hair, as hair may become entangled in the rotating massage heads.
- Use on injured, irritated, or swollen skin, or on areas with open wounds.
- Use for medical purposes, including treatment of diseases or injuries.
- Use by children or persons who are unable to safely operate the device.
- Continuous use on the same body area for longer than recommended.
- Use of the device while charging.
- Use of the device under running water/immersion of the main unit in water – **please refer to the safety information for clear instructions.**
- Use of the device if it is damaged or malfunctioning.
- Modification, disassembly, or repair of the device by unauthorized persons.
- Looking directly into the LED lights during operation.

The manufacturer is not liable for damage caused by improper use.

APPLICATION

STEP 1: PREPARATION

Ensure the device is fully charged before use. The power button flashes red during charging, then turns green once fully charged.

Select the desired massage head. Attach the massage head securely to the device by aligning the middle and pressing down. Always ensure that the device is switched off before attaching or removing accessories.

To remove the massage head, grab it firmly and lift it off the connecting pin.

NOTE: If you do not use the device for a longer period, we recommend charging the device at least every three months to avoid damage to the battery.


STEP 2: OPERATION

Press the **power button** to turn the device on.

Select the **desired massage mode**:

- Forward rotation with red LED lights
- Backward rotation with red LED lights
- Forward rotation
- Backward rotation

Adjust the speed using the “strength” button according to your preference.

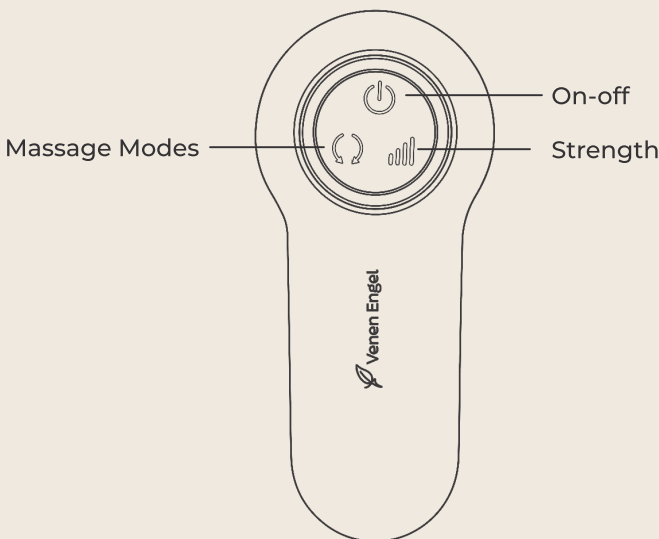Apply the massage head gently to the desired body area and move the device slowly across the skin. Take your time to discover your personal comfort level. The application should always feel soothing and must never be uncomfortable or painful. Please ensure that individual sessions do not exceed 15 minutes per body area, and no longer than 30 minutes of continuous use overall.

Avoid directing the LED light toward the eyes or looking straight into the light.

STEP 3: AFTER USE

Once your session has finished, **press the power button** to turn off the device. To remove a massage head, grab it firmly and lift it off the connecting pin.

CLEANING & STORAGE

CLEANING

To keep your Venen Engel 3-in-1 Glow Brush in perfect condition, we recommend cleaning it regularly. Always switch off the device before cleaning. Clean the main unit using a soft, slightly damp cloth. Do not use aggressive cleaning agents. The removable massage heads may be cleaned separately under running water if necessary. Ensure that all parts are completely dry before reassembling, and do not allow water to enter the main unit.

STORAGE AND TRANSPORTATION

We recommend storing your device in a dry, clean, and dust-free environment. Avoid direct sunlight, excessive heat, and high humidity. Transport the device in its original packaging whenever possible.

SAFETY INSTRUCTIONS

Keep this instruction manual in a safe place! Read these instructions for use carefully before using the device and use it only as described in this manual. The device is intended for capable adults only; keep the device and its accessories out of the reach of children. Do not use the device if it or the USB charging cable is damaged or malfunctioning. Do not modify, disassemble, or repair the device yourself, as repairs may only be carried out by authorized specialists, otherwise the warranty will be void.

The device is intended solely for cosmetic and wellness purposes and is not a medical device. Do not use it for medical purposes, including the treatment of diseases or injuries. Consult a physician before use if you are pregnant, have a serious heart condition, or have injured, irritated, or swollen skin, or open wounds in the treatment areas. Stop using the device immediately if you feel discomfort, pain, or notice any unusual skin reactions.

Do not use the device on the head or near hair to avoid entanglement, pulling, or pinching. Do not look directly into the integrated red LED lights during operation, and do not direct the light into the eyes of other persons. The LED light is intended solely as part of the device's function and is not intended for therapeutic or medical use.

To ensure product longevity and safety the main unit should not be exposed to direct water jets, continuous running water, or immersed in water. Remove the massage heads to clean them separately with water, and ensure all parts are completely dry before reattaching. Do not use the device while it is charging, and only charge it using the supplied USB cable. Always turn off the device before changing the massage heads. Do not operate the device continuously for more than 30 minutes, and do not use it on the same body area for more than 15 minutes to avoid skin sensitivity issues.

WARRANTY

We offer a 24-month warranty on the control unit. In the event of malfunction, it may be necessary to send the device back for inspection. Please ensure that the shipping box is padded to avoid possible damage during transport. Unfortunately, no claims under the warranty can be brought for defects caused by improper shipping. Insofar as a legal warranty obligation exists, a replacement will either be provided free of charge or the manufacturer shall have the option to repair the device, under exclusion of the right of cancellation or a reduction in the purchase price. If repeated attempts to effect repair are unsuccessful, or a replacement unit also proves to be defective due to the manufacturer's fault, the customer shall be entitled to cancel the purchase or receive a reduction in the purchase price.

Warranty exclusion: Excluded from the warranty are defects due to damage caused by the use of force, improper operation, externally applied force, or modifications and repair work performed by third parties, e.g. defects caused by incorrectly rated or short-circuited fuses, or defects attributable to normal wear and tear.

Warranty requirements: The warranty claim exists only in connection with the original proof of purchase (invoice). Please keep the proof of purchase in a safe place.

Should you have any problems or questions when using the system, please do not hesitate to contact us by mail at service@venenengel.de.

DISPOSAL INSTRUCTIONS

IMPORTANT: Please do not dispose of the product in your household waste!

GERMANY:

Electronic products must not be disposed of in unsorted municipal waste. Instead, use the available return and collection systems. In this way, we can all contribute to the reuse, recycling, and recovery of electronic products.

Incorrect disposal of electronic products not only leads to a waste of resources (most electrical appliances contain valuable raw materials such as gold, copper, silver, etc.), but also endangers the environment due to the sometimes hazardous substances they contain.

ALTERNATIVE: We are happy to take care of the disposal of your old Venen Engel device. Simply send it to:

**Wellcosan GmbH
Bellinger Tor 16a
36396 Steinau a. d. Str.
Germany**

Please enclose a note requesting disposal.

Old batteries and accumulators must be disposed of separately from solid municipal waste. You can do this using the return and recycling systems available to you (e.g., Corepile).

Improper disposal of old batteries and accumulators can have a significant impact on humans and the environment due to the hazardous substances they contain. Recycling preserves valuable raw materials and allows product components to be reused.

EUROPEAN UNION:

Electronic products must not be disposed of in unsorted municipal waste. Instead, use the available return and collection systems. This way, we can all contribute to the reuse, recycling, and recovery of electronic products.

Incorrect disposal of electronic products not only wastes resources (most electrical appliances contain valuable raw materials such as gold, copper, silver, etc.), but also endangers the environment due to the sometimes hazardous substances they contain.

Used batteries and accumulators must be disposed of separately from solid municipal waste. You can do this using the return and recycling systems available to you (e.g., Corepile).

Improper disposal of old batteries and accumulators can have a significant impact on people and the environment due to the hazardous substances they contain. Recycling preserves valuable raw materials and allows product components to be reused.

LIABILITY NOTICE

If you have health concerns or are unsure, please consult a healthcare professional before use.

NOTE:

Your Venen Engel product is **not a medical device** and is not intended for self-treatment.

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Wellcosan GmbH, c/o Caya Postbox 794627,