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OPERATING INSTRUCTIONS

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
OPERATING INSTRUCTIONS

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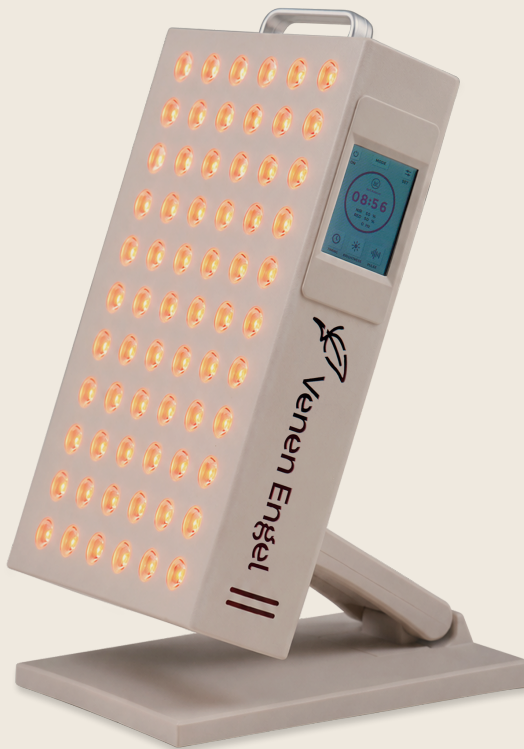PRELIMINARY REMARK

Congratulations on purchasing this device! You will enjoy using your Venen Engel InfraGlow Product! If you are satisfied with your Venen Engel, we would be delighted if you would recommend it to others.

This manual contains details that users should observe in order to prevent hazards and use the device safely.

If you have any questions about use, accessories, or suggestions, please do not hesitate to contact us!

Mail service@venenengel.de

Web www.venenengel.de

BEFORE YOU START

- Please read the operating instructions carefully before using the device.
- Keep this instruction manual in a safe place so that you can refer to it when necessary.
- The Venen Engel product must not be used in certain health conditions. Please read the list of contraindications carefully.

IMPORTANT SAFETY INSTRUCTIONS AND WARNINGS

It is important that you read all warnings and precautions in this manual. They are for your safety, to prevent injury, and to avoid damage to the device.

- Please read these instructions for use carefully before using the device.
- Keep this manual for future reference.
- Use the device only as described in this manual.
- Do not use the device if it is damaged.
- The device is intended for adults who are able to understand the instructions. Persons with reduced physical, sensory, or mental capabilities should only use the device under supervision.
- Keep the device away from children.
- The device is intended for indoor use only in dry environments. Do not use the device in bathrooms or other humid environments.
- Do not use the device for medical purposes or therapeutic claims.
- Never use the light overnight.
- Wear eye-protection goggles at all times during application. Do not look directly into the light or the LEDs.
- Do not use the device if you feel discomfort during use. Stop using the device immediately if pain or unusual discomfort occurs.
- Please consult a physician before using the device if you are pregnant.
- Consult a physician before use if you have a known sensitivity to light or are taking photosensitizing medication.
- Consult a physician before use if you suffer from conditions that may be triggered by light exposure (e.g., certain forms of epilepsy).
- Consult a physician before use if you have injured, irritated, inflamed, or infected skin.
- Consult a physician before use if you have reduced heat sensitivity or cannot properly perceive heat.
- Do not use the device if the power cable is damaged. Ensure you operate the device only with the specified voltage (AC 100-240V).

- Do not operate the device with wet hands.
- Do not immerse the device in water or other liquids.
- Disconnect the device from the power supply before cleaning.
- Do not cover the device during operation.
- A session of 10–30 minutes, 3–5 times a week at a distance of approx. 20–50 cm is recommended. Do not exceed 30 minutes per session.
- Allow the panel to cool down before storing it.
- Do not modify, disassemble, or repair the device. The panel does not contain user-serviceable parts, and repairs must be carried out by qualified personnel only.

SCOPE OF DELIVERY

EW 300 light panel with stand

Eye protection goggles

Remote control

Connecting cable

TECHNICAL DATA

Modell	EW300
Wavelengths	630 nm 660 nm 810 nm 850 nm
LEDs	60pcs 5W Quad-Chip
Input	AC 100-240 V, 50/60 HZ
Power consumption	approx. 110 W
Dimensions	12.3 x 9 x 2.6 inches (31 x 23 x 6.6 cm)
Weight	3.3kg
Protection Class	Class I
IP Rating	IP20
Brightness Range	10 - 100%
Time Setting	1 - 90 minutes
Operating Temperature	10 - 60°C
Operating Humidity	< 85 %
Storage Temperature	10 - 60°C


EXPLANATION OF SYMBOLS

	Manufacturer information
	Made in China
	Conforms with EU standards
	Use only indoors
	Electronic product. Dispose according to instructions
	Read manual before use
	Serial number
	Protect from rain
	Outer packaging made of carton
	Inner packaging made of polypropylene
	Disposal instructions France for electronic products
	Disposal instructions France for the packaging

PRODUCT DESCRIPTION

This product is designed to support general wellness, relaxation, and daily self-care routines through the use of red and near-infrared light exposure. Specially engineered for precise, targeted applications, this compact light panel brings an intensive wellness ritual right into your home as a non-medical device. Its versatile design allows the panel to be mounted, making it easy to seamlessly integrate into your favorite cozy corner.

The product utilizes a harmonized combination of visible red light (630nm & 660nm) and near-infrared light (810nm & 850nm). While the soothing red light is ideal for a dedicated moment of mindfulness, leaving your skin and senses feeling revitalized and cared for, the near-infrared light—which is completely invisible to the human eye—promotes gentle relaxation and comfort. It is the perfect ritual to help your mind and body wind down and de-accelerate after a long day.

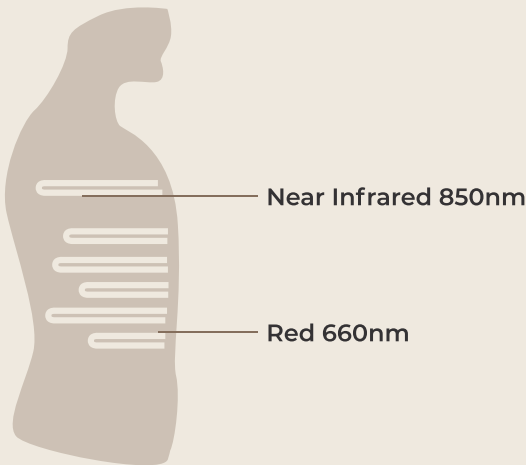INTENDED USE

The **Venen Engel InfraGlow Mini** is intended for general wellness through exposure to red and near-infrared light. It may support general well-being, relaxation, and a pleasant feeling in the body during and after use. The device is designed for use on the body surface in a non-invasive manner.

This product is **not a medical device** and is not intended for diagnosis, treatment, or prevention of any disease.

TARGET USER

The device is intended for:

- Adults
- Individuals able to understand the instructions

Persons with reduced physical, sensory or mental capabilities should only use the device under supervision.

CONTRAINDICATIONS

PLEASE CONSULT A PHYSICIAN BEFORE USING THE DEVICE IF YOU:

- are pregnant
- have known sensitivity to light or are taking photosensitizing medication
- have injured, irritated, inflamed or infected skin
- suffer from conditions that may be triggered by light exposure (e.g. certain forms of epilepsy)
- have reduced heat sensitivity or cannot properly perceive heat.

If you are unsure if the device is suitable for you, consult a physician before use. Stop using the device immediately if you feel discomfort or pain.

IMPROPER USE

Any use other than the intended use is considered improper. Examples of foreseeable misuse include:

- Use for medical purposes or therapeutic claims
- Use on injured, infected, or irritated skin
- Direct eye exposure without protection
- Covering the device during operation
- Use in wet or humid environments (e.g. bathroom)
- Exceeding the recommended application time

The manufacturer is not liable for damage caused by improper use.

APPLICATION

STEP 1: BEFORE USE

Place the panel securely on a flat, stable surface and adjust its stand to choose the desired angle. Ensure that the panel is installed securely. Connect the cable to a power outlet. Flip the green power switch on the back of the panel. Put on the safety goggles (or keep them within reach) and turn your back towards the panel, so the lights do not hit your eyes turning the device on.

STEP 2: OPERATION

- **Flip the power switch at the back of the panel to turn it on.**
- **Either press the power button on the remote control or press the power button on the touch control panel on the side of the panel to switch it completely on.**
- **Select the desired mode via the remote control or by using the touch panel at the side of the panel.** You can select one of the eight pre-built programs or create a custom program.
- **Adjust the brightness** of each frequency (between 10% - 100%) within the program. You can increase or decrease the brightness in 10% steps.
- **Adjust the pulse frequency (between 0 - 100 Hz).** 0Hz means no pulse, 100Hz means very quick pulse. You can increase or decrease the frequency in 1Hz-steps.
- **Set the application time (1-90 minutes).** We recommend a session of 10-30 minutes, 3-5 times a week at a distance of approx. 20-50 cm. Please ensure that individual sessions should not exceed 30 minutes.

Important: Always wear the supplied protective eyewear during every use, regardless of the selected settings or program. Never look directly into the LEDs. This also applies when using near-infrared light, even though it is invisible to the human eye. Regular sunglasses do not provide adequate protection and must not be used as a substitute for the supplied protective eyewear.

PROGRAM / GUIDE

With the **Venen Engel InfraGlow Mini**, you can treat yourself to conscious wellness breaks exactly where you need them most. This compact LED light panel proves that true power is not a matter of size. Specially designed for precise, targeted application, it brings deep relaxation, better inner balance, and a soothing moment of mindfulness to your daily self-care routine. It is the perfect, flexible companion to help you easily integrate a targeted light ritual into your busy day.

AREAS OF APPLICATION

The handily sized InfraGlow Mini is designed for focused, point-by-point application on specific body areas. Whether placed on your desk, nightstand, or in your favorite cozy corner, it serves as a personal oasis of calm to help you switch off and practice conscious breathing.

- **Soothed skin and senses:** The comforting visible red light (630 & 660 nm) provides a pampered body feel, giving your skin a conscious moment of mindfulness and leaving it with a refreshed, natural radiance.
- **Cozy relaxation:** The invisible, gentle near-infrared light (810 & 850 nm) offers deep, soothing comfort, helping you to ease your mind and de-accelerate after a long, demanding day.

DURATION AND FREQUENCY

The recommended session time is **10 to 20 minutes** (up to a maximum of 30 minutes), ideally **3 to 5 times a week**. Take your time to discover your personal comfort level: you can fully customize your experience by adjusting the brightness (**10–100%**) and the pulse frequency (**0–100 Hz**), or by selecting one of our pre-set wellness programs like Glow Boost, Deep Comfort, or your own tailored My Routine starting from 10 minutes.

Remember, more is not always better—the consistency of your new light ritual is the absolute key to achieving a fantastic, well-balanced body feeling. The application should always feel soothing and must never be uncomfortable or painful. Longer sessions do not necessarily improve results. Consistency is the secret to the best experience.

PROGRAMS

VE	NIR	RED	Time in min
Skin Radiance	0%	100%	20
Vitalize & Activate	100%	0%	15
Deep Comfort	100%	0%	30
Renew & Recharge	0%	100%	25
Glow Boost	100%	100%	15
Soft Balance	50%	50%	15
My Routine	customize your own programm		from 10
Soothe & Reset	100%	0%	20
Evening Wind-down	0%	50%	30

1. SKIN RADIANCE (20 MIN.)

Light: 100% RED

This 20-minute beauty ritual for your skin and senses invites you to enjoy a moment of mindfulness and treat your skin to an all-round nourished body feeling. The red light gently caresses your skin surface, softly bringing out your natural glow. Red wavelengths in the range of 630–660 nm are among the most thoroughly researched spectra in dermatological science.

2. VITALIZE & ACTIVATE (15 MIN.)

Light: 100% NIR

This 15-minute vitalizing ritual gently awakens your body and provides you with fresh energy for the day. The invisible near-infrared light (NIR) envelops your body with a gentle glow and boosts your energetic body feeling. From a scientific perspective, near-infrared light in the spectrum of 810–850 nm penetrates particularly deep into the tissue.

3. DEEP COMFORT (30 MIN.)

Light: 100% NIR

This 30-minute wellness ritual offers you an intense timeout of pure recovery and invites you to experience deep relaxation. Just like Vitalize & Activate, Deep Comfort relies exclusively on 100% near-infrared light (NIR). Due to the longer application time of the Deep Comfort program, the light gently caresses your body and harmonizes your overall well-being.

4. RENEW & RECHARGE (25 MIN.)

Light: 100% RED

This 25-minute regeneration ritual invites you to recharge your batteries and give your skin an intensive self-care break. It focuses 100% on powerful red light that caresses your skin surface, bringing your skin tone into a fresh, harmonious balance. Compared to Skin Radiance, the longer application duration here ensures an even deeper sense of relaxation.

5. GLOW BOOST (15 MIN.)

Light: 100% NIR | 100% RED

This 15-minute intensive ritual offers a highly efficient upgrade to your daily beauty routine, providing instant freshness in no time. It combines invisible near-infrared light with powerful red light. This ideal combination makes your complexion truly glow.

6. SOFT BALANCE (15 MIN.)

Light: 50% NIR | 50% RED

This 15-minute harmony ritual invites you to pause during a stressful day, clear your senses, and find a new balance. The secret of this program lies in its perfectly balanced light composition: both wavelength ranges – near-infrared and red light – work together gently and harmoniously. This combination softly caresses your senses and brings your mind to a state of relaxed calm.

7. MY ROUTINE (AB 10 MIN.)

Light: Individually customizable

This highly flexible program invites you to design your me-time entirely freely and tailor it precisely to your form on any given day. You decide exactly how to mix the proportions of near-infrared light (NIR) and red light to suit your needs.

8. SOOTHE & RESET (20 MIN.)

Light: 100% NIR

This 20-minute relaxation ritual invites you to take a deep breath after a demanding day and gently let excess tension melt away. The program relies entirely on the pure power of invisible near-infrared light. It envelops your body in deep-acting, invigorating energy and brings your entire sense of physical well-being into a state of deep rest.

9. EVENING WIND-DOWN (30 MIN.)

Light: 50% RED

This 30-minute evening ritual invites you to gently leave the day's stimuli behind and carefully prepare your mind for a peaceful night's rest. The program relies on a softly dosed, pure red light at 50% intensity. This reduced intensity gently caresses your skin in the evening hours, smoothly guiding your entire system into relaxation mode.

AFTER USE

TERMINATION AND DISCONNECTION

Once your session has finished, switch off the device using either the remote control or the touch control panel, then turn off the green power switch on the back of the device. You can also stop the session at any time by pressing the power button.

After use, we recommend disconnecting the panel from the power supply immediately when it is not in use. Before storing the panel, please make sure to let it cool down completely.

CLEANING & STORAGE

CLEANING

To keep your InfraGlow Mini in perfect condition, we recommend cleaning it regularly. Always make sure to disconnect the device from the power supply before you begin cleaning. We recommend using a dry or slightly damp cloth to gently wipe the surface area of the panel. Please do not use any aggressive cleaning agents, chemical solvents, or harsh detergents, as they may damage the material.

MAINTENANCE

Please note that the InfraGlow Mini does not contain any user-serviceable parts. To ensure your safety and prevent any damage, repairs must be carried out by qualified personnel only.

STORAGE AND TRANSPORTATION

We recommend storing your InfraGlow Mini in a dry, clean, and dust-free environment. Make sure to avoid exposing the device to direct sunlight and strictly adhere to the recommended storage temperature range of 10°C to 60°C.

TRANSPORTATION

If you need to transport your InfraGlow Mini, we recommend using the original packaging whenever possible. This ensures the device is optimally protected from physical impact and shock during transit.

SAFETY INSTRUCTIONS

Keep this operating manual in a safe place! Read these instructions for use carefully before using the device and use it only as described in this manual. Do not use the device if it is damaged; if damage is detected, unplug it and switch it off immediately. Do not modify, dismantle, or repair the device yourself. Repairs to the device may only be carried out by specialists authorized by Wellcosan GmbH, otherwise the warranty will be void. **Wear eye-protection goggles at all times during application and do not look directly into the light.** If you feel discomfort, or if pain or unusual discomfort occurs, stop using the device immediately. The device may become warm during use. It must not be covered by other objects during operation, nor should it be used on heat-sensitive surfaces.

Only use the device with the specified power supply and ensure that the power supply voltage matches the device specifications. Do not operate the device or handle the power plug with wet hands. Do not use the device in bathrooms or other humid environments, and never immerse the device in water or other liquids. If the device has been exposed to or come into contact with liquids, disconnect it immediately from the power supply. When not in use, always unplug the power cord. Always unplug the device from the power outlet before cleaning. Keep the product away from children or persons who have not been instructed in its use.

WARRANTY

We offer a 24-month warranty on the device. In the event of malfunction, it may be necessary to send the device back for inspection. Please ensure that the shipping box is padded to avoid possible damage during transport. Unfortunately, no claims under the warranty can be brought for defects caused by improper shipping. Insofar as a legal warranty obligation exists, a replacement will either be provided free of charge or the manufacturer shall have the option to repair the device, under exclusion of the right of cancellation or a reduction in the purchase price. If repeated attempts to effect repair are unsuccessful, or a replacement unit also proves to be defective due to the manufacturer's fault, the customer shall be entitled to cancel the purchase or receive a reduction in the purchase price.

Warranty exclusion: Excluded from the warranty are defects due to damage caused by the use of force, improper operation, externally applied force, or modifications and repair work performed by third parties, e.g. defects caused by incorrectly rated or short-circuited fuses, or defects attributable to normal wear and tear.

Warranty requirements: The warranty claim exists only in connection with the original proof of purchase (invoice). Please keep the proof of purchase in a safe place.

Should you have any problems or questions when using the system, please do not hesitate to contact us by mail at service@venenengel.de.

DISPOSAL INSTRUCTIONS

IMPORTANT: Please do not dispose of the product in your household waste!

GERMANY

Electronic products must not be disposed of in unsorted municipal waste. Instead, use the available return and collection systems. In this way, we can all contribute to the reuse, recycling, and recovery of electronic products.

Incorrect disposal of electronic products not only wastes resources (most electrical appliances contain valuable raw materials such as gold, copper, silver, etc.), but also endangers the environment due to the sometimes hazardous substances they contain.

Alternative: We are happy to take care of the disposal of your old Venen Engel device. Simply send it to:

Wellcosan GmbH
Bellinger Tor 16a
36396 Steinau a. d. Str.
Germany

Please enclose a note requesting disposal.

EUROPEAN UNION

Electronic products must not be disposed of in unsorted municipal waste. Instead, use the return and collection systems available. In this way, we can all contribute to the reuse, recycling, and recovery of electronic products.

Incorrect disposal of electronic products not only wastes resources (most electrical appliances contain valuable raw materials such as gold, copper, silver, etc.), but also endangers the environment due to the sometimes hazardous substances they contain.

LIABILITY NOTICE

If you have health concerns or are unsure, please consult a healthcare professional before use.

NOTE

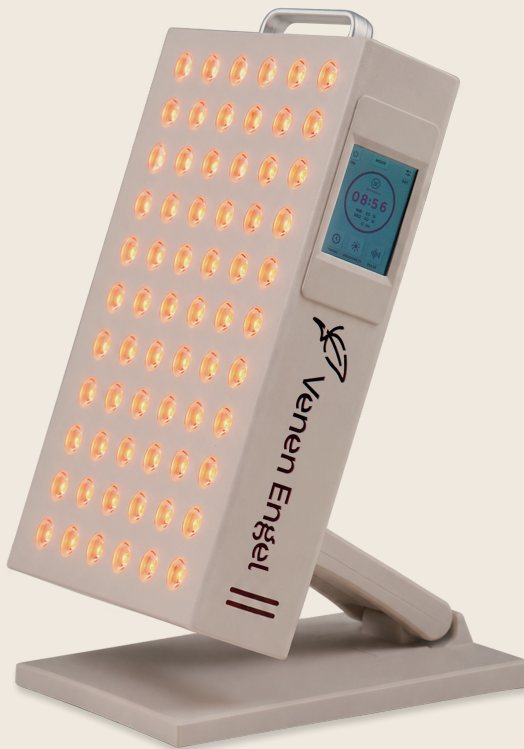
Your Venen Engel product is **not a medical device** and is not intended for self-treatment.

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© Wellcosan GmbH 2026
Wellcosan GmbH, c/o Caya Postbox 794627,