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OPERATING INSTRUCTIONS

Venen Engel InfraGlow Pro

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Venen Engel InfraGlow Pro

PRELIMINARY REMARK

Congratulations on purchasing this device! You will enjoy using your Venen Engel InfraGlow Product! If you are satisfied with your Venen Engel, we would be delighted if you would recommend it to others.

This manual contains details that users should observe in order to prevent hazards and use the device safely.

If you have any questions about use, accessories, or suggestions, please do not hesitate to contact us!

Mail service@venenengel.de

Web www.venenengel.de

BEFORE YOU START

- Please read the operating instructions carefully before using the device.
- Keep this instruction manual in a safe place so that you can refer to it when necessary.
- The Venen Engel product must not be used in certain health conditions. Please read the list of contraindications carefully.

IMPORTANT SAFETY INSTRUCTIONS AND WARNINGS

It is important that you read all warnings and precautions in this manual. They are for your safety, to prevent injury, and to avoid damage to the device.


- Please read these instructions for use carefully before using the device.
- Keep this manual for future reference.
- Use the device only as described in this manual.
- Do not use the device if it is damaged.
- The device is intended for adults who are able to understand the instructions. Persons with reduced physical, sensory, or mental capabilities should only use the device under supervision.
- Keep the device away from children.
- The device is intended for indoor use only in dry environments. Do not use the device in bathrooms or other humid environments.
- Do not use the device for medical purposes or therapeutic claims.
- Never use the light overnight.
- Wear eye-protection goggles at all times during application. Do not look directly into the light or the LEDs.
- Do not use the device if you feel discomfort during use. Stop using the device immediately if pain or unusual discomfort occurs.
- Please consult a physician before using the device if you are pregnant.
- Consult a physician before use if you have a known sensitivity to light or are taking photosensitizing medication.
- Consult a physician before use if you suffer from conditions that may be triggered by light exposure (e.g., certain forms of epilepsy).
- Consult a physician before use if you have injured, irritated, inflamed, or infected skin.
- Consult a physician before use if you have reduced heat sensitivity or cannot properly perceive heat.

- Do not use the device if the power cable is damaged. Ensure you operate the device only with the specified voltage (AC 100-240V).
- Do not operate the device with wet hands.
- Do not immerse the device in water or other liquids.
- Disconnect the device from the power supply before cleaning.
- Do not cover the device during operation.
- A session of 10–30 minutes, 3–5 times a week at a distance of approx. 20–50 cm is recommended. Do not exceed 30 minutes per session.
- Allow the panel to cool down before storing it.
- Do not modify, disassemble, or repair the device. The panel does not contain user-serviceable parts, and repairs must be carried out by qualified personnel only.

SCOPE OF DELIVERY

EW1000 light panel

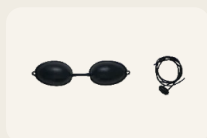Eye protection goggles

Remote control

Connecting cable

Cable and snap links

Door hook

Hanging hardware

TECHNICAL DATA

Model	EW1000
Wavelengths	480nm 525nm 630nm; 660nm 810nm 830nm; 850nm 940nm 980nm; 1060nm
LEDs	200 pcs 5W Dual Chip
Input	AC 100-240 V, 50/60 Hz
Power consumption	approx. 307 W
Dimensions	36 x 9.1 x 2.6in (91 x 23 x 6.6 cm)
Weight	6.52 kg
Protection class	Class I
IP rating	IP20
Brightness range	10 – 100%
Time setting	1 – 90 minutes
Operating temperature	10 – 60 °C
Operating humidity	< 85%
Storage temperature	10 – 60 °C


EXPLANATION OF SYMBOLS

	Manufacturer information
	Made in China
	Conforms with EU standards
	Use only indoors
	Electronic product. Dispose according to instructions
	Read manual before use
	Serial number
	Protect from rain
	Outer packaging made of carton
	Inner packaging made of Polypropylene
	Disposal instructions France for electronic products
	Disposal instructions France for the packaging

PRODUCT DESCRIPTION

This product is designed to support general wellness, relaxation, and daily self-care routines through a multi-dimensional light experience. Specially engineered for precise, targeted applications, this light panel brings an intensive wellness ritual right into your home as a non-medical device. Featuring **8 preset wellness programs** and a customizable **"My Routine"** mode, its versatile design allows the panel to be mounted with the included door hook or used with an optional, rollable stand—making it easy to seamlessly integrate into your favorite cozy corner.

The product utilizes a harmonized combination of visible **red light (630nm & 660nm)**, near-infrared light (**810nm to 1060nm**), **blue light (480nm)**, and **green light (525nm)**, which dynamically change in intensity depending on your chosen program. While the soothing red light is ideal for a dedicated moment of mindfulness, leaving your skin and senses feeling revitalized and cared for, the invisible near-infrared light delivers gentle warmth to promote deep relaxation and a pleasant body sensation. Complemented by blue light for refreshing clarity and green light for a gentle balance, it is the perfect ritual to help your mind and body wind down and de-accelerate after a long day.

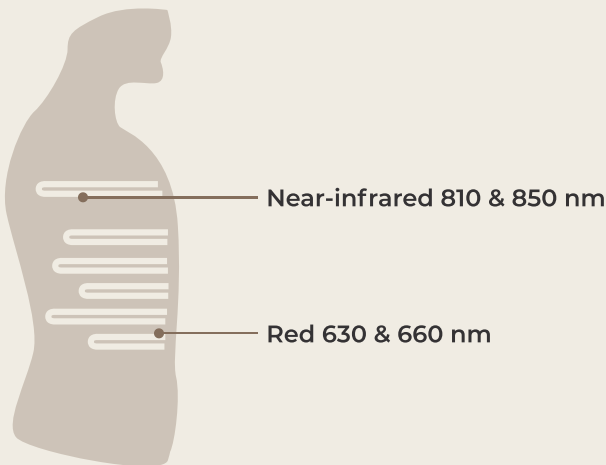INTENDED USE

The **Venen Engel InfraGlow Pro** is intended for general wellness through exposure to red and near-infrared light. It may support general well-being, relaxation, and a pleasant feeling in the body during and after use. The device is designed for use on the body surface in a non-invasive manner.

This product is **not a medical device** and is not intended for diagnosis, treatment, or prevention of any disease.

CONTRAINDICATIONS

PLEASE CONSULT A PHYSICIAN BEFORE USE IF YOU:

- are pregnant
- have known sensitivity to light or are taking photosensitizing medication
- have injured, irritated, inflamed or infected skin
- suffer from conditions that may be triggered by light exposure (e.g. certain forms of epilepsy)
- have reduced heat sensitivity or cannot properly perceive heat

If you are unsure if the device is suitable for you, consult a physician before use. **Stop using the device immediately if you feel discomfort or pain.**

TARGET USER

The device is intended for:

- Adults
- Individuals able to understand the instructions

Persons with reduced physical, sensory or mental capabilities should only use the device under supervision.

IMPROPER USE

Any use other than the intended use is considered improper.

Examples of foreseeable misuse include:

- Use for medical purposes or therapeutic claims
- Use on injured, infected, or irritated skin
- Direct eye exposure without protection
- Covering the device during operation
- Use in wet or humid environments (e.g. bathroom)
- Exceeding the recommended application time


The manufacturer is not liable for damage caused by improper use.

APPLICATION


ASSEMBLING YOUR PRODUCT

If you want to hang your InfraGlow Pro to your wall or door, what you need:


- ① Threaded posts
- ② Cable support
- ③ Door hook
- ④ Cable gripper
- ⑤ Snap link

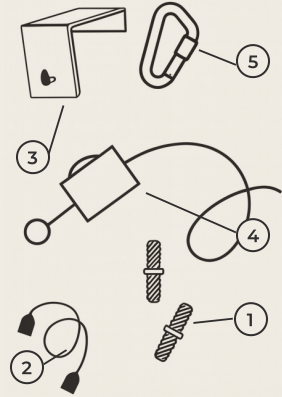
- Screw the threaded posts into the top of your InfraGlow Pro device.
- Screw the cable supports to the top of the threaded posts.

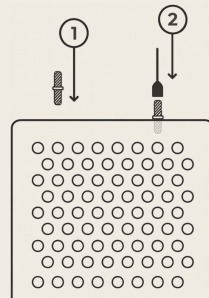- Hang the door hook on the top of any standard-sized door.

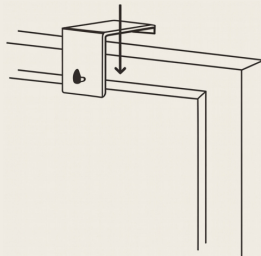- Connect the cable supports with the snap-link.
- Hang the other end of the cable directly onto the door hook. The snap link and attached cable gripper allows you to adjust the height of your panel individually.

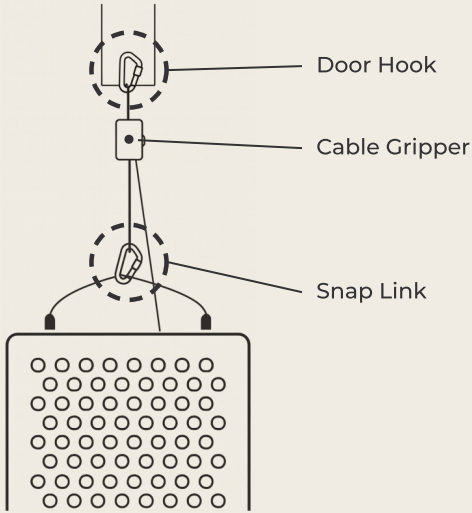IF YOU WANT TO MOUNT YOUR RED LIGHT PANEL ON THE STAND

- Move the base to where you would like the panel to stand.
- Carefully slot the panel between the two columns of the stand.

STEP 1: BEFORE USE

Set up the panel flexibly to suit your needs: hang it securely using the included door hook, place it on a stable surface, preferably leaning slightly against a wall, or use the optional rolling stand—always ensuring the device is safely installed and maintains the recommended minimum distance.

Connect the cable to a power outlet and switch on the power on the back of the panel. For optimal safety and maximum relaxation, always put on the included safety goggles and keep your eyes closed throughout the entire application as you comfortably start your ritual using the intuitive control panel or remote control.

STEP 2: OPERATION

- **Flip the power switch** at the back of the panel to turn it on.
- **Either press the power button on the remote control or press the power button on the touch control panel** on the side of the panel to switch it completely on.
- **Select the desired mode via the remote control or by using the touch panel** at the side of the panel. You can select one of the eight pre-built programs or create a custom program.
- **Adjust the brightness** of each frequency (between 10% – 100%) within the program. You can increase or decrease the brightness in 10% steps.
- **Adjust the pulse frequency (between 0 – 100 Hz)**. 0Hz means no pulse, 100Hz means very quick pulse. You can increase or decrease the frequency in 1Hz steps.
- **Set the application time (1–90 minutes)**. We recommend a session of 10–30 minutes, 3–5 times a week at a distance of approx. 20–50 cm. Please ensure that individual sessions do not exceed 30 minutes.

PROGRAM / GUIDE

With the **Venen Engel InfraGlow Pro**, you can treat yourself to conscious wellness breaks exactly where you need them most. Designed for large-area applications, this versatile LED light panel brings deep relaxation, a gentle balance, and a soothing moment of mindfulness to your daily self-care routine. It is the perfect, flexible companion to help you easily integrate a tailored light ritual for your skin and senses into your busy day.

AREAS OF APPLICATION

The flexibly adaptable InfraGlow Pro is designed for a completely customized wellness experience. Whether hung safely using the included door/wall hook, set up in your favorite cozy corner, or moved freely through the room on the optional rollable stand, it serves as a personal oasis of calm to help you switch off and enjoy a moment of pure mindfulness.

- **Soothed skin and senses:** the comforting visible red light (630 & 660 nm) provides a pampered body feel, giving your skin a conscious moment of mindfulness and leaving it feeling refreshed, natural and radiant.
- **Refreshing and balancing:** the InfraGlow Pro also features blue light (480 nm) for refreshing clarity and green light (525 nm) for a gentle balance.
- **Cozy relaxation:** the invisible, gentle near-infrared light (810 & 850 nm) offers deep, soothing comfort, helping you to ease your mind and de-accelerate after a long, demanding day.

DURATION AND FREQUENCY

The recommended session time is 10 to 20 minutes (up to a maximum of 30 minutes), ideally 3 to 5 times a week. Take your time to discover your personal comfort level: you can fully customize your experience by adjusting the brightness (10–100%) and the pulse frequency (0–100 Hz), or by selecting one of our pre-set wellness programs like Glow Boost, Deep Comfort, or your own tailored My Routine starting from 10 minutes.

Remember, more is not always better—the consistency of your new light ritual is the absolute key to achieving a fantastic, well-balanced body feeling. The application should always feel soothing and must never be uncomfortable or painful. Longer sessions do not necessarily improve results. Consistency is the secret to the best experience.

PROGRAMS

	NIR	RED	BLUE	GREEN	Time in min
Skin Radiance	0%	100%	50%	50%	20
Vitalize & Activate	100%	0%	0%	0%	15
Deep Comfort	100%	0%	0%	0%	30
Renew & Recharge	0%	100%	100%	0%	25
Glow Boost	100%	100%	100%	0%	15
Soft Balance	50%	50%	50%	50%	15
My Routine	customize your own program				from 10
Soothe & Reset	100%	0%	100%	0%	20
Evening Wind-down	0%	50%	0%	50%	30

1. Skin Radiance (20 Min.)

Light: 100% RED | 50% BLUE | 50% GREEN | 0% NIR

This 20-minute beauty ritual for your skin and senses invites you to enjoy a moment of mindfulness and treat your skin to an all-around nourished feeling. Powerful red light is accompanied by gentle hints of blue and green light. This combination caresses the surface of your skin, softly illuminating your natural radiance. Red wavelengths in the 630–660 nm range are among the most extensively researched spectra in dermatological science.

2. Vitalize & Activate (15 Min.)

Light: 100% NIR | 0% RED | 0% BLUE | 0% GREEN

This 15-minute vitalizing ritual gently awakens your body and provides you with renewed energy for the day. The invisible, near-infrared light (NIR) caresses your body with a gentle glow, invigorating your physical energy and well-being. From a scientific perspective, near-infrared light in the 810–850 nm spectrum penetrates particularly deep into the tissue.

3. Deep Comfort (30 Min.)

Light: 100% NIR | 0% RED | 0% BLUE | 0% GREEN

This 30-minute well-being ritual treats you to an intensive timeout of pure relaxation and invites you to experience deep tranquility. Just like Vitalize & Activate, Deep Comfort relies exclusively on 100% near-infrared light (NIR). Thanks to the longer application time of the Deep Comfort program, the light caresses your body and brings your entire well-being into perfect harmony.

4. Renew & Recharge (25 Min.)

Light: 100% RED | 100% BLUE | 0% NIR | 0% GREEN

This 25-minute regeneration ritual invites you to recharge your batteries and treat your skin to a vitalizing skincare break. The program relies on a dual light composition: powerful red light (100%) and clarifying blue light (100%) work hand in hand. This combination caresses the surface of your skin and brings your complexion into a fresh, harmonious balance.

5. Glow Boost (15 Min.)

Light: 100% NIR | 100% RED | 100% BLUE | 0% GREEN

This 15-minute intensive ritual treats you to a highly efficient upgrade for your daily beauty routine and instantly lends you a newfound freshness. It combines invisible near-infrared light, powerful red light, and clarifying blue light into an intensive triad. This combination makes your complexion glow.

6. Soft Balance (15 Min.)

Light: 50% NIR | 50% RED | 50% BLUE | 50% GREEN

This flexible program invites you to shape your me-time with absolute freedom and tailor it precisely to how you feel each day. You decide exactly how to mix the proportions of near-infrared light (NIR), red light, blue light, and green light.

7. My Routine (from 10 Min.)

Light: individually customizable

This highly flexible program invites you to design your me-time entirely freely and tailor it precisely to your form on any given day. You decide exactly how to mix the proportions of near-infrared light (NIR) and red light to suit your needs.

8. Soothe & Reset (20 Min.)

Light: 100% NIR | 0% RED | 100% BLUE | 0% GREEN

This 20-minute relaxation ritual invites you to take a deep breath after a demanding day and gently let excess tension melt away. The program pairs the invisible, near-infrared light directly with clarifying blue light. This combination caresses your body with invigorating light and surface freshness, bringing your entire physical well-being into a state of deep rest.

9. Evening Wind-down (30 Min.)

Light: 0% NIR | 50% RED | 0% BLUE | 50% GREEN

This 30-minute evening ritual invites you to wind down the day in an oasis of tranquility and gently prepare your body for a restful night. A balanced mix of 50% red and 50% green light envelops you in a gentle, soothing colored-light bath. This combination gently brings your inner rhythm into harmony.

AFTER USE

TERMINATION AND DISCONNECTION

Once your session has finished, please switch off the device both with the remote control or via the touch control panel and by the green switch on the back of the product. You can also stop the application at any time by pressing the power button.

After use, we recommend disconnecting the panel from the power supply immediately when it is not in use. Before storing the panel, please make sure to let it cool down completely.

CLEANING & STORAGE

CLEANING

To keep your Venen Engel InfraGlow Pro in perfect condition, we recommend cleaning it regularly. Always make sure to disconnect the mat from the power supply before you begin cleaning. We recommend using a dry or slightly damp cloth to gently wipe the surface area of the panel. Please do not use any aggressive cleaning agents, chemical solvents, or harsh detergents, as they may damage the material.

MAINTENANCE

Please note that the mat does not contain any user-serviceable parts. To ensure your safety and prevent any damage, repairs must be carried out by qualified personnel only.

STORAGE AND TRANSPORTATION

We recommend storing your Venen Engel InfraGlow Pro in a dry, clean, and dust-free environment. Make sure to avoid exposing the device to direct sunlight and strictly adhere to the recommended storage temperature range of 10–60°C.

TRANSPORTATION

If you need to transport your InfraGlow Pro, we recommend using the original packaging whenever possible. This ensures the device is optimally protected from physical impact and shock, during transit.

SAFETY INSTRUCTIONS

Keep this operating manual in a safe place! Read these instructions for use carefully before using the device and use it only as described in this manual. Do not use the device if it is damaged; if damage is detected, unplug it and switch it off immediately. Do not modify, dismantle, or repair the device yourself. Repairs to the device may only be carried out by specialists authorized by Wellcosan GmbH, otherwise the warranty will be void. Wear eye-protection goggles at all times during application and do not look directly into the light. If you feel discomfort, or if pain or unusual discomfort occurs, stop using the device immediately. The device may become warm during use. It must not be covered by other objects during operation, nor should it be used on heat-sensitive surfaces.

Only use the device with the specified power supply and ensure that the power supply voltage matches the device specifications. Do not operate the device or handle the power plug with wet hands. Do not use the device in bathrooms or other humid environments, and never immerse the device in water or other liquids. If the device has been exposed to or come into contact with liquids, disconnect it immediately from the power supply. When not in use, always unplug the power cord. Always unplug the device from the power outlet before cleaning. Keep the product away from children or persons who have not been instructed in its use.

WARRANTY

We offer a 24-month warranty on the control unit. In the event of malfunction, it may be necessary to send the device back for inspection. Please ensure that the shipping box is padded to avoid possible damage during transport. Unfortunately, no claims under the warranty can be brought for defects caused by improper shipping. Insofar as a legal warranty obligation exists, a replacement will either be provided free of charge or the manufacturer shall have the option to repair the device, under exclusion of the right of cancellation or a reduction in the purchase price. If repeated attempts to effect repair are unsuccessful, or a replacement unit also proves to be defective due to the manufacturer's fault, the customer shall be entitled to cancel the purchase or receive a reduction in the purchase price.

Warranty exclusion: Excluded from the warranty are defects due to damage caused by the use of force, improper operation, externally applied force, or modifications and repair work performed by third parties, e.g. defects caused by incorrectly rated or short-circuited fuses, or defects attributable to normal wear and tear.

Warranty requirements: The warranty claim exists only in connection with the original proof of purchase (invoice). Please keep the proof of purchase in a safe place.

Should you have any problems or questions when using the system, please do not hesitate to contact us at +49 (0) 666 174 891 05 or by mail at service@venenengel.de.

DISPOSAL INSTRUCTIONS

IMPORTANT: Please do not dispose of the product in your household waste!

GERMANY

Electronic products must not be disposed of in unsorted municipal waste. Instead, use the available return and collection systems. In this way, we can all contribute to the reuse, recycling, and recovery of electronic products.

Incorrect disposal of electronic products not only wastes resources (most electrical appliances contain valuable raw materials such as gold, copper, silver, etc.), but also endangers the environment due to the sometimes hazardous substances they contain.

Alternative: We are happy to take care of the disposal of your old Venen Engel device. Simply send it to:

Wellcosan GmbH
Bellinger Tor 16a
36396 Steinau a. d. Str.
Germany

Please enclose a note requesting disposal.

OTHER COUNTRIES / EU

Electronic products must not be disposed of in unsorted municipal waste. Instead, use the return and collection systems available. In this way, we can all contribute to the reuse, recycling, and recovery of electronic products.

Incorrect disposal of electronic products does not only waste resources (most electrical appliances contain valuable raw materials such as gold, copper, silver, etc.), but also endangers the environment due to the sometimes hazardous substances they contain.

LIABILITY NOTICE

If you have health concerns or are unsure, please consult a healthcare professional before use.

NOTE:

Your Venen Engel product is **not a medical device** and is not intended for self-treatment.

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© Wellcosan GmbH 2026
Wellcosan GmbH, c/o Caya Postbox 794627,